

Task Number #156641

Testing ordered > 07 MAY '26

Sample received > 12 MAY '26

Client drakonpeptides.com

Sample Selank 10mg

Manufacturer Unknown

Batch DRK-SEL-MAY26

Sample description >

See picture or pictures attached.

Tests requested >

Selank analysis

Results >



Selank	9.28 mg
Purity	99.616%

Comments >

Sample was tested twice.

Analysis conducted > 14 MAY 2026

Signature > 





SHIMADZU

LabSolutions


Analysis Report

<Sample Information>

Sample Name : 20948-1 Selank 10mg, LCMS Screening
 Sample ID : 052
 Data Filename : 20948-1 Selank 10mg, LCMS Screening_060.lcd
 Method Filename : 251007_ACN-H2O_Peptidy2_0.5.lcm
 Batch Filename : 20260513_2 - PO VYPNUTI ELEKTRINY.lcb
 Vial # : 4-50 Sample Type : Unknown
 Injection Volume : 0.2 uL
 Date Acquired : 14.05.2026 1:50:35 Acquired by : System Administrator
 Date Processed : 14.05.2026 12:21:44 Processed by : System Administrator

<Chromatogram>

mAU

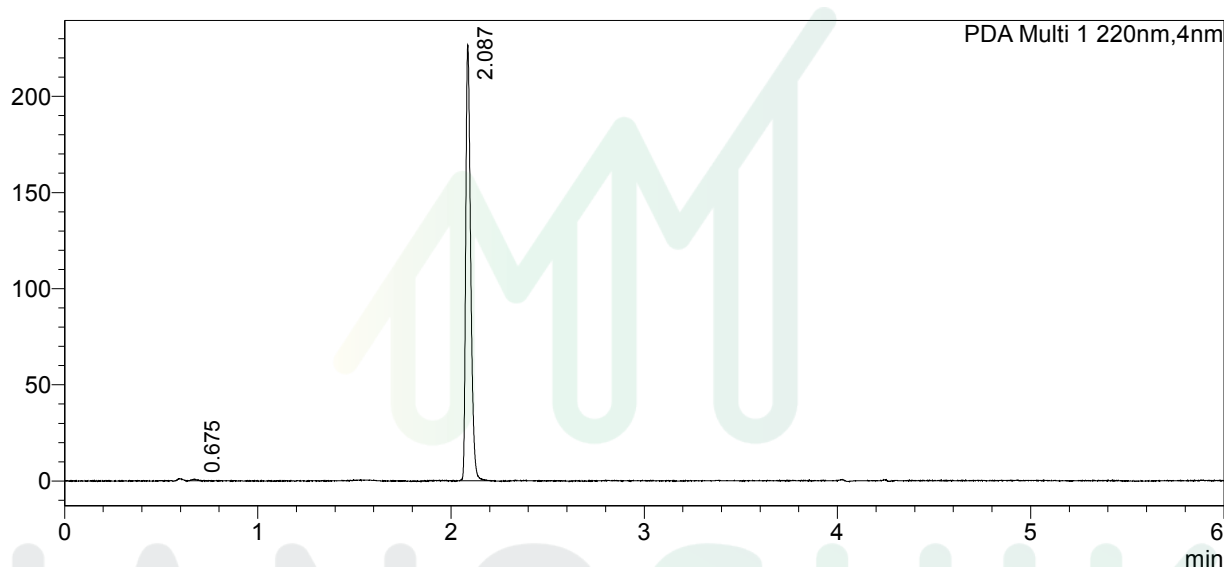

<Peak Table>

PDA Ch1 220nm

Peak#	Ret. Time	Area	Conc.	Unit	Name	Area%
1	0.675	1558	0.00			0.384
2	2.087	404500	9.28	mg	Selank	99.616
Total		406059				100.000

Task Number #156643

Testing ordered > 07 MAY '26

Sample received > 12 MAY '26

Client drakonpeptides.com

Sample Selank 10mg

Manufacturer Unknown

Batch DRK-SEL-MAY26

Sample description >

See picture or pictures attached.

Tests requested >

Heavy Metals Analysis

Results >



As	not detected.
Cd	not detected.
Pb	not detected.

Comments >

Hg: not detected. Although the raw data show low ppm values for heavy metals, these do not necessarily indicate true detection. As shown in the attached graphs, the signals are consistent with electronic noise inherent to the detection system; therefore, the results were concluded as “not detected”.

Analysis conducted > 20 MAY 2026

Signature > 

Verify this test at www.janoshik.com/verify/ with the following unique key

AQ7S8WANE4TS




Sample : 20948-1
 Operator: MB
 Comment : Comment
 Group : As_test
 Date : 2026-05-20 10:56:37

Measurement Condition

Instrument	Atmosphere	Collimator	Sample Cup
EDX-7000	Air	10 (mm)	Polypropylene

Analyte	TG kV	uA	FI	Acq. (keV)	Anal. (keV)	Time (sec)	DT (%)
As	Rh 50	17-Auto	----	0 - 40	10.30-10.80	Live- 50	29

[cps/uA]

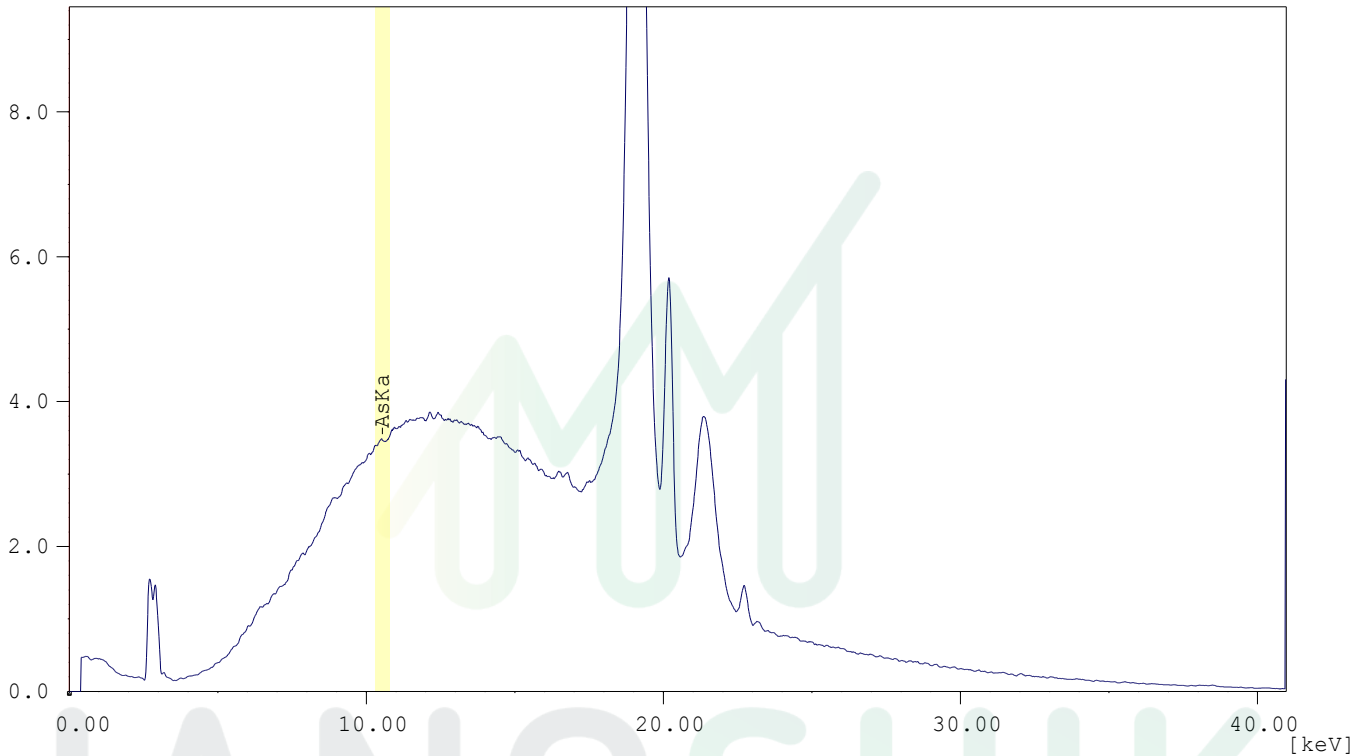

Peak List

Channel	Line	keV	Net Int. (cps/uA)
As	AsKa	10.54	0.1600

Quantitative Result

Analyte	Result	[3-sigma]	Proc.-Calc.	Line	Int. (cps/uA)
As	0.7 ppm	L [4.1]	Quan-EC	AsKa	0.1600

Sample : 20948-1
 Operator: MB
 Comment : Comment
 Group : Hg_Pb_Cd
 Date : 2026-05-20 10:34:07

Measurement Condition

Instrument: EDX-7000 Atmosphere: Air Collimator: 10 (mm) Sample Cup: Polypropylene

Analyte	TG kV	uA	FI	Acq. (keV)	Anal. (keV)	Time (sec)	DT (%)
Cd	Rh 50	640-Auto	#5	0 - 40	22.72-23.52	Live- 50	30
Hg	Rh 50	640-Auto	#5	0 - 40	9.74-10.24	Live- 50	30
Pb	Rh 50	640-Auto	#5	0 - 40	12.38-12.88	Live- 50	30

Qualitative Result

Element: Cd, Hg, Pb

Peak List

Channel	Line	keV	Net Int. (cps/uA)
Cd	CdKa	23.12	0.0007 Q
Hg	HgLa	9.98	0.0000 Q
Pb	PbLb1	12.62	0.0017 Q

Quantitative Result

Analyte	Result	[3-sigma]	Proc.-Calc.	Line	Int. (cps/uA)
Cd	3.9 ppm	L [3.5]	Quan-EC	CdKa	0.0007
Hg	0.2 ppm	L [1.2]	Quan-EC	HgLa	0.0000
Pb	3.4 ppm	L [4.5]	Quan-EC	PbLb1	0.0017

JANOSHNIK


Sample : 20948-1
Operator: MB

Group : Hg_Pb_Cd
Comment : Comment


Date : 2026-05-20 10:31:07
01/01




[cps/uA] Cd

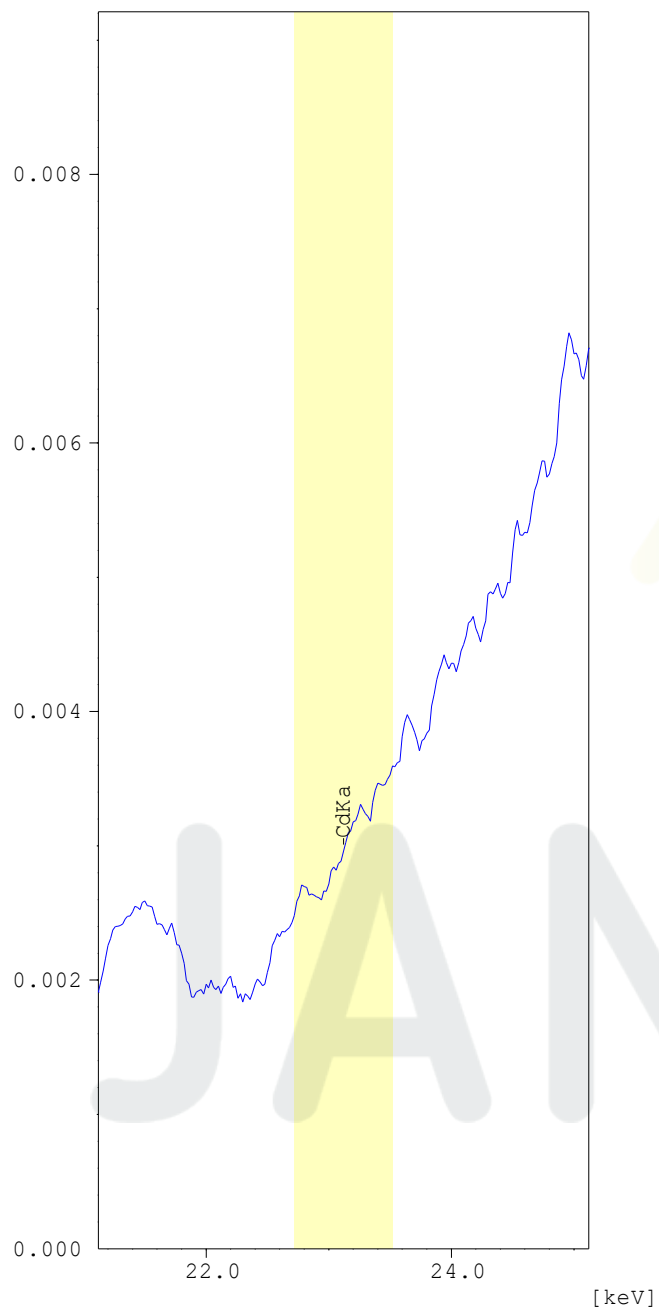

[cps/uA] Hg

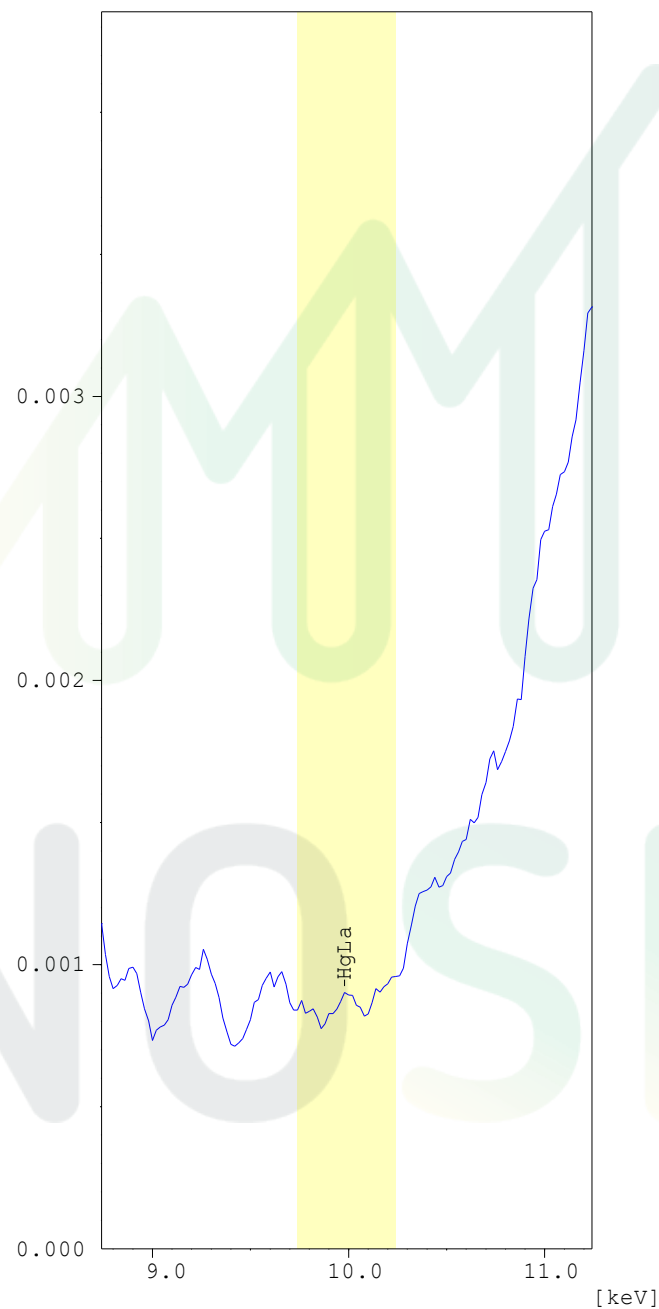

[cps/uA] Pb

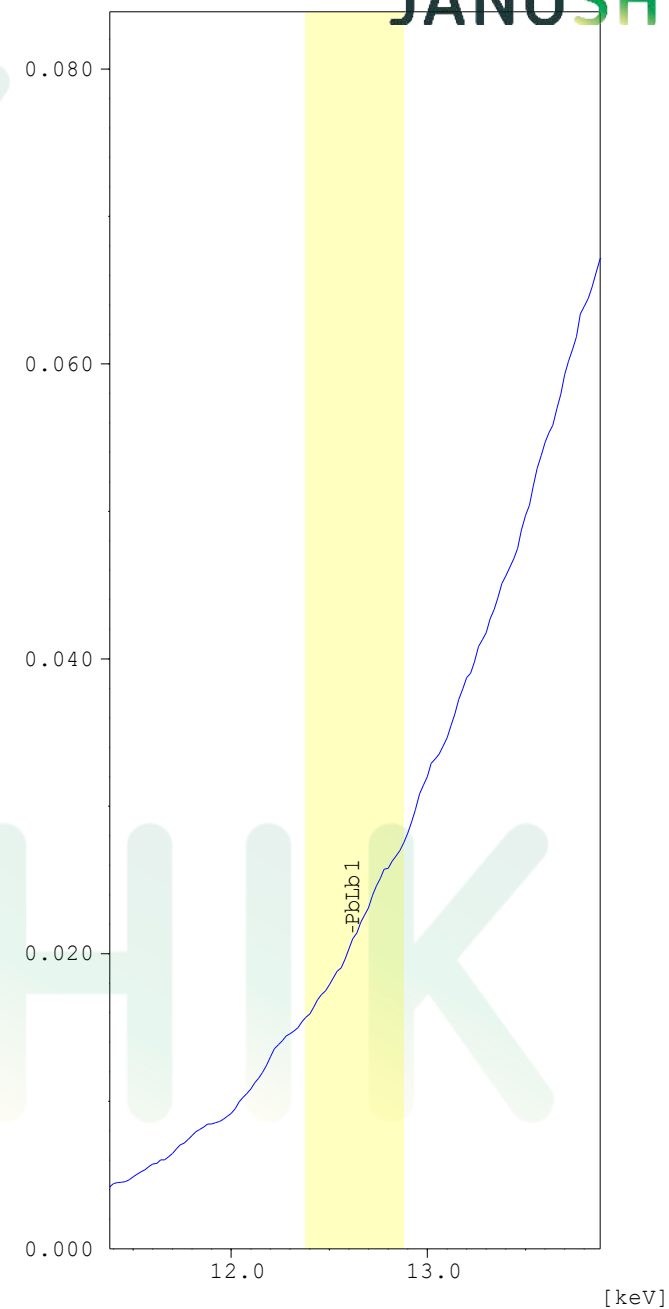

Task Number #156642

Testing ordered > 07 MAY '26

Sample received > 12 MAY '26

Clientdrakonpeptides.com

SampleSelank 10mg

ManufacturerUnknown

BatchDRK-SEL-MAY26

Sample description >

See picture or pictures attached.

Tests requested >

LCMS Screening For Peptide Contamination

Results >



LCMS screening	

Comments >

LCMS analysis identified the presence of mannitol, selank, and degradation products. LCMS spectra are attached. No significant contamination detected.

Analysis conducted > 18 MAY 2026

Signature > 

Verify this test at www.janoshik.com/verify/ with the following unique key

M9I5KW2A8PCA






JANOSHIK

==== Shimadzu LabSolutions Data Report ====

<Spectrum>

Line#:1 R.Time:0.264(Scan#:124)
MassPeaks:73
RawMode:Single 0.264(124) BasePeak:205(7354681)
BG Mode:None Segment 1 - Event 1

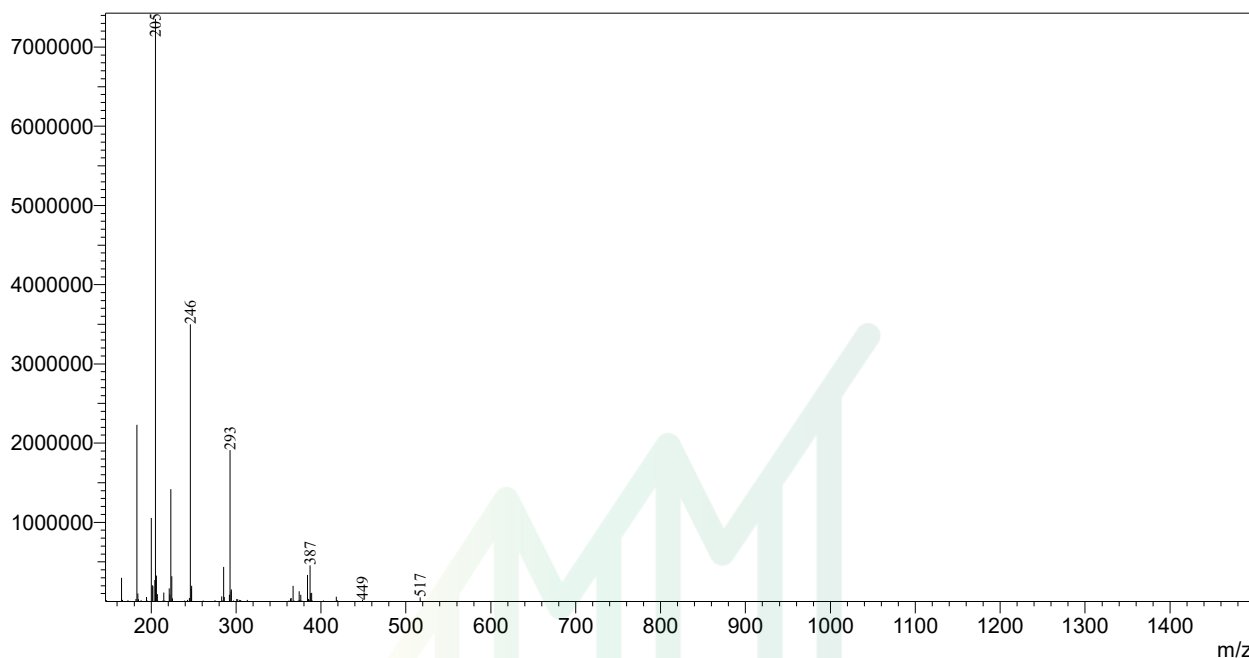




JANOSHIK

==== Shimadzu LabSolutions Data Report ====

<Spectrum>

Line#:1 R.Time:0.357(Scan#:167)
MassPeaks:116
RawMode:Single 0.357(167) BasePeak:377(16734656)
BG Mode:None Segment 1 - Event 1

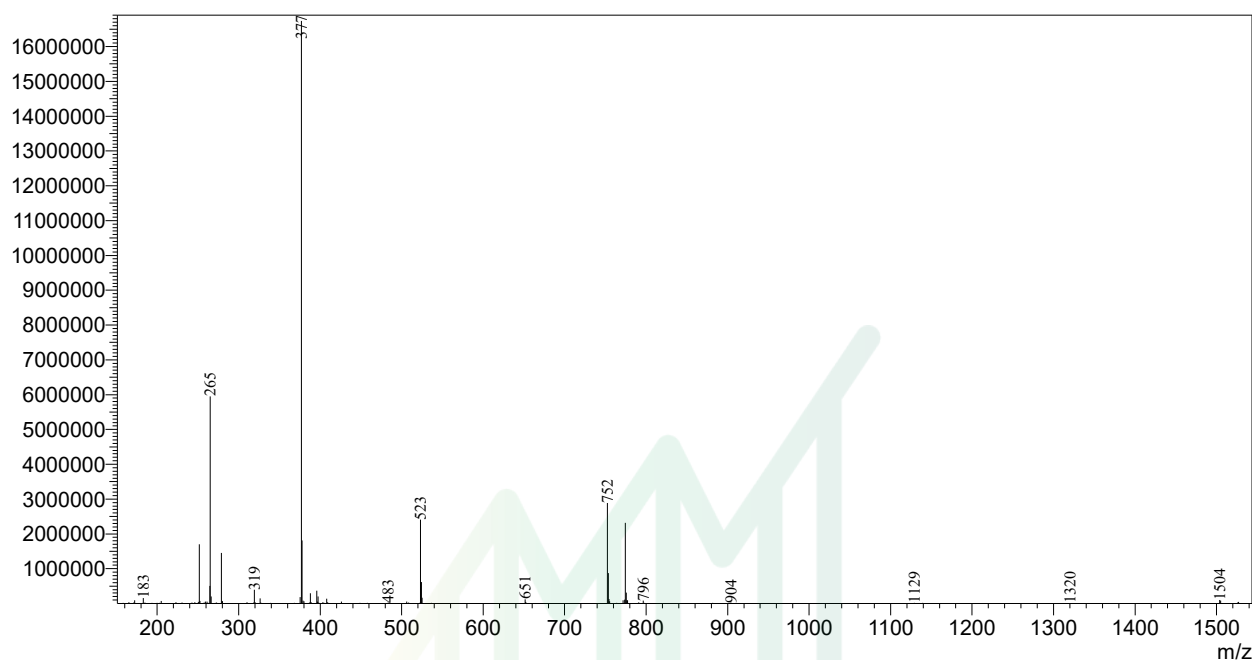




JANOSHIK

==== Shimadzu LabSolutions Data Report ====

<Spectrum>

Line#:1 R.Time:0.387(Scan#:181)
MassPeaks:95
RawMode:Single 0.387(181) BasePeak:377(13920911)
BG Mode:None Segment 1 - Event 1

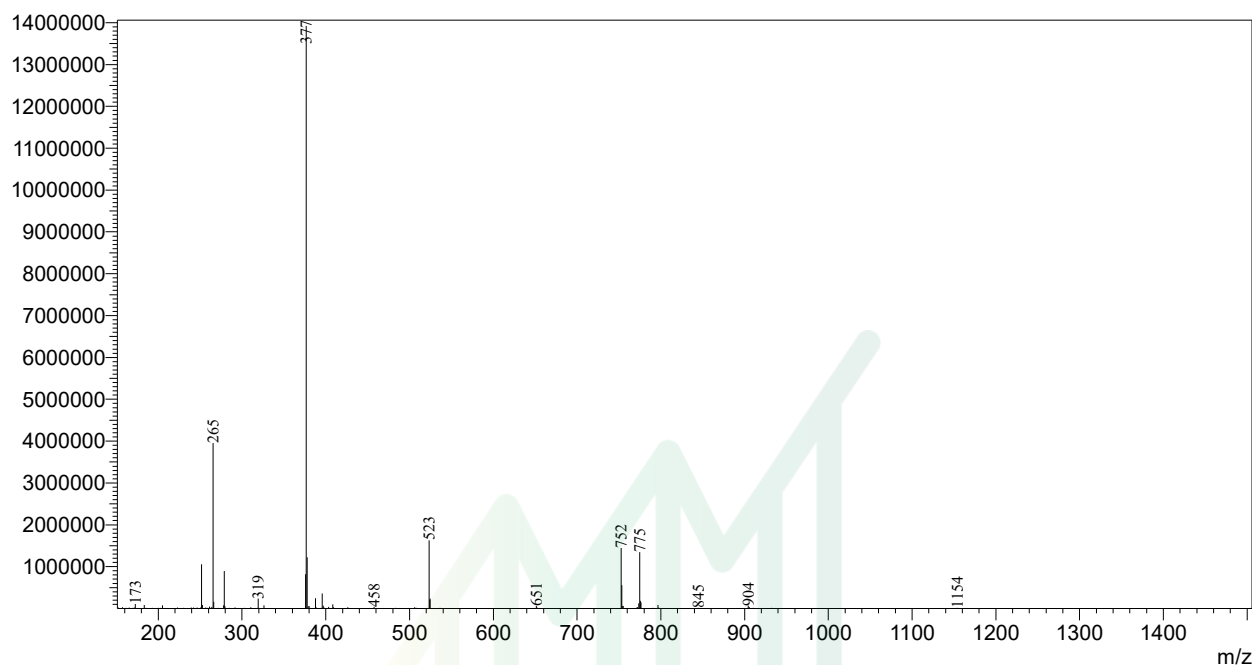




JANOSHIK

==== Shimadzu LabSolutions Data Report ====

<Chromatogram>

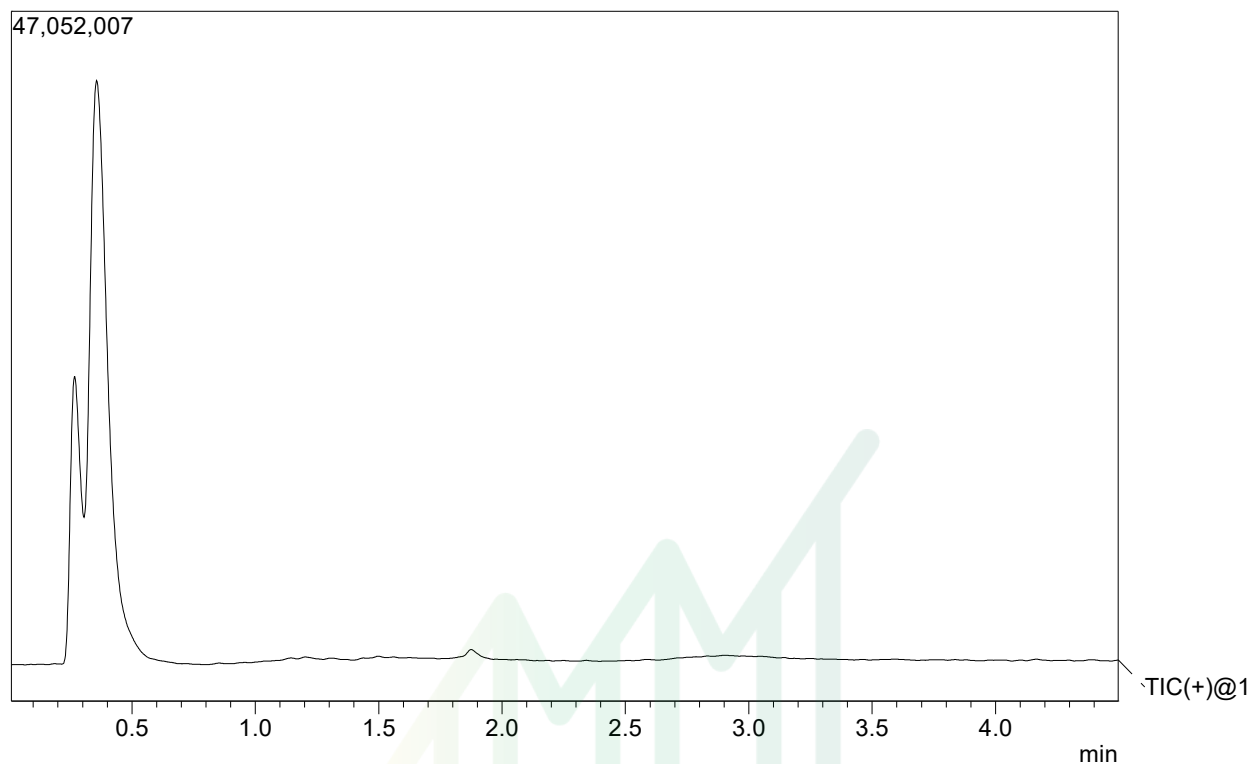

JANOSHIK